


Tracto Gastrointestinal

Endoprótesis Autoexpandible para Tracto Gastrointestinal

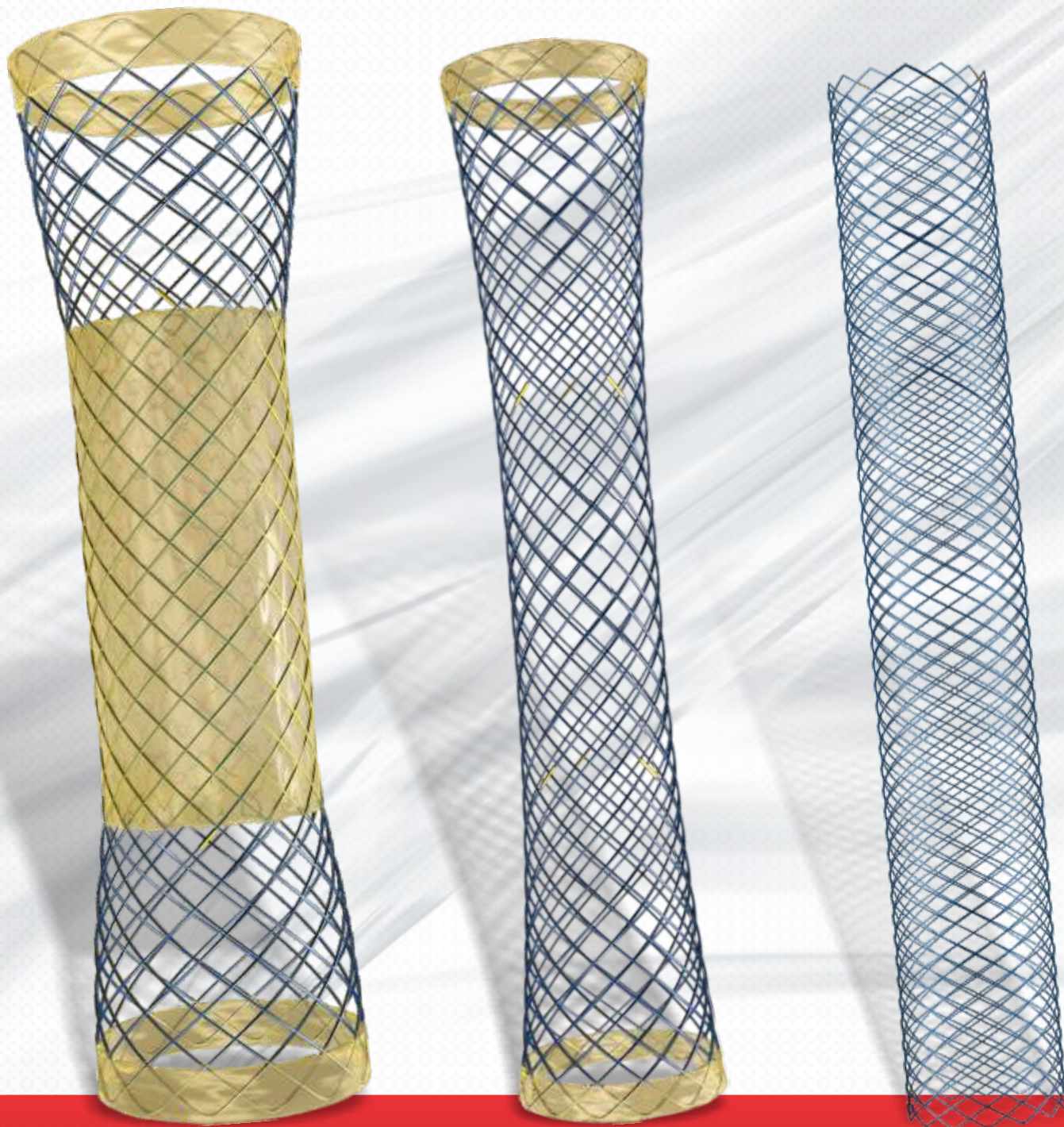

www.braile.com.br

 **BRAILE**[®]
BIOMÉDICA

Modelos

Biliar



08 x 50
08 x 75
08 x 100
10 x 50
10 x 75
10 x 100

Biliar revestida



08 x 50
08 x 75
10 x 50
10 x 75

Esofágico não revestido



20 x 70
20 x 100
20 x 130

Colon



26 x 70
26 x 100
26 x 130

Esofágica revestida



20 x 70
20 x 100
20 x 130

Enteral



20 x 70
20 x 100
20 x 130

Catéteres


Biliares
08 Fr | 1900mm

Biliares Revestidos
10 Fr | 1900mm

Enterais
10 Fr | 2500mm

Colon y Esofágicas
23 Fr | 700mm

La Endoprótesis Autoexpansible para Tracto Gastrointestinal es constituida de hilos de Nitinol, aleación metálica superelástica y memoria térmica, tranzada por un único hilo (monofilamento), en forma tubular formando rombos (diamond shape) y tiene la borda distal acabada con marcadores de oro. La prótesis posee revestimiento parcial de poliuretano.

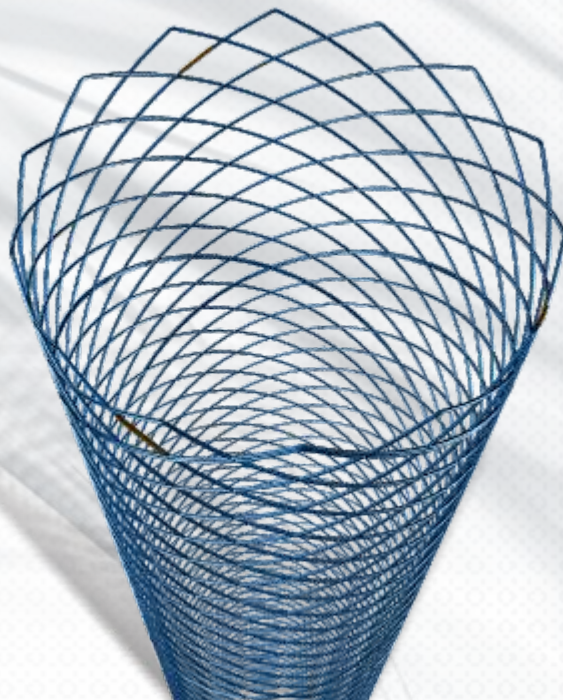

Braile Biomédica Indústria, Comércio e Representações Ltda.

Av. Pres. Juscelino K. de Oliveira, 1505 • Jd. Tarraf I
Código Postal 15091-450 • São José do Rio Preto • SP • Brasil
Teléfono (55) 17 2136 7000 • Fax 17 2136 7040
Other countries • Phone # (55) 17 2136 7000 • Fax 17 2136 7040
CNPJ: 52.828.936/0001-09
Técnico Responsable: Vladimir D. A. Ramirez • CRF/SP 9010
Registro ANVISA nº 10159030069
www.braille.com.br

BRAILE
BIOMÉDICA