


Organic Tubular Arterial Graft Valvulated Tube 300 Smooth

Experience

Quality

Performance

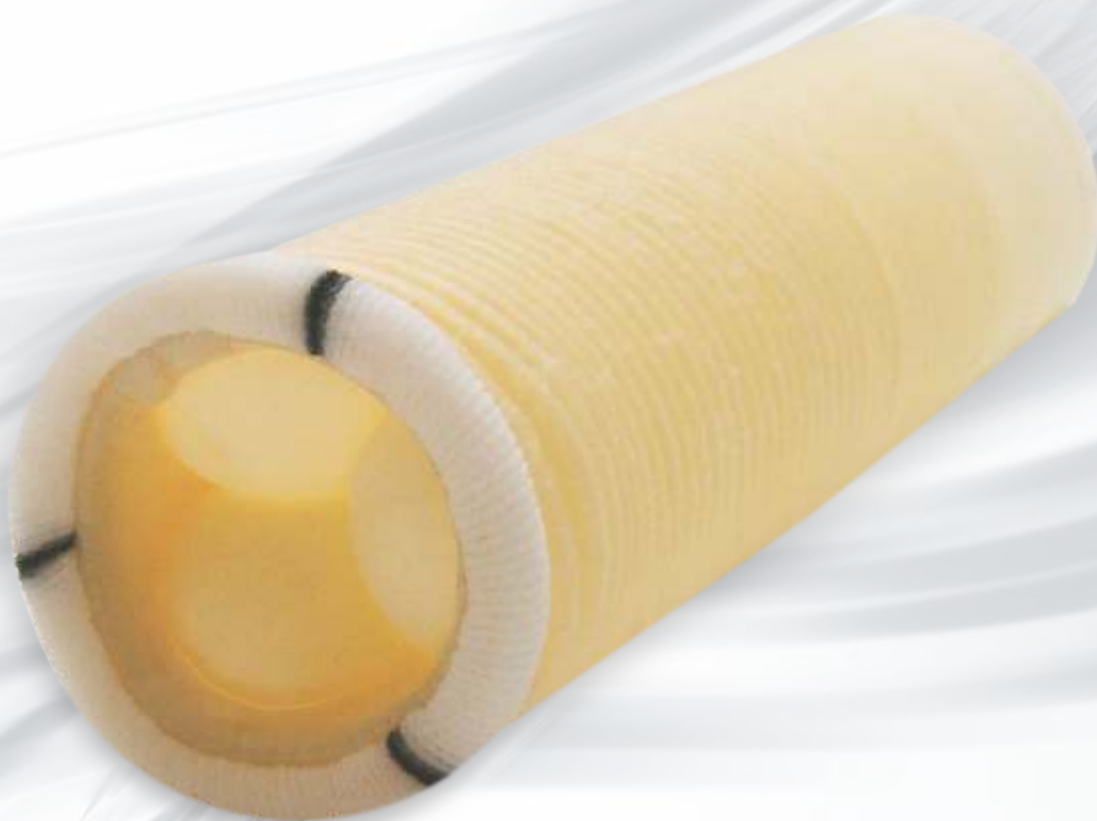


www.braile.com.br

 **BRAILE**[®]
BIOMÉDICA

Description: The Organic Tubular Arterial Graft Valvulated Tube 300 Smooth is made of bovine pericardium treated with glutaraldehyde and conserved in formaldehyde 4%. This treatment gives it proper characteristics of resistance, flexibility and lack of antigenicity.

Description: Organic Tubular Arterial Graft Valvulated Tube 300 Smooth is indicated for corrections of aneurysms of the ascending aorta with valvar insufficiency

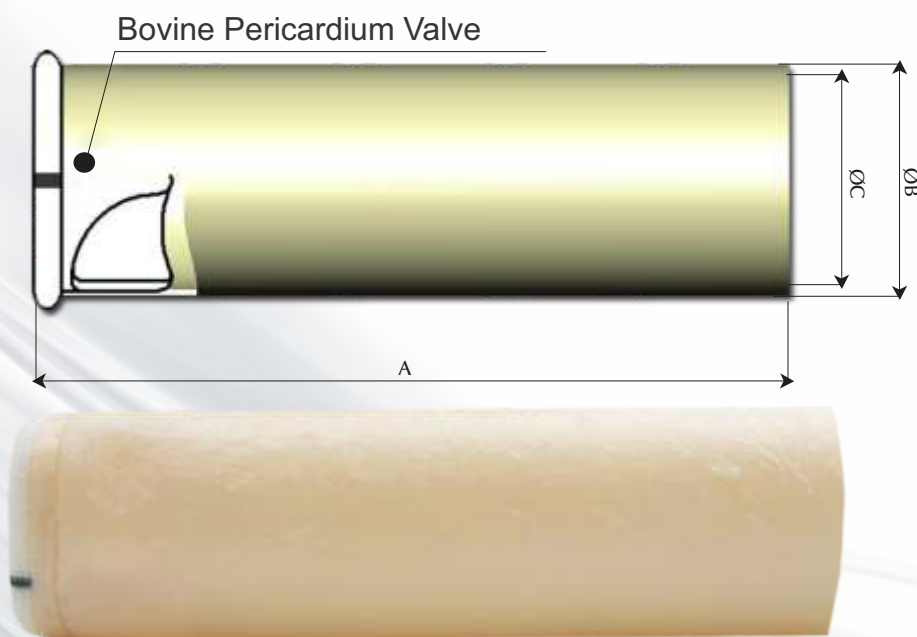

Technical Specifications (mm)

Size	A ± 10	ØB ± 0,5	ØC ± 0,5	Valve	Code
25	120	27	25	Ø25	10935
27	120	29	27	Ø27	10942
29	120	31	29	Ø29	10959
31	120	33	31	Ø31	10966
33	120	35	33	Ø33	10973

The size of the tube is specified by the valve Ø.

Performance Information

Diameter	Average Pressure Grad. (mmHg)	Effective Orifice Area (cm ²)
25	4,73	2,38
27	3,00	2,55
29	1,53	3,08
31	1,32	3,63
33	0,56	3,88



Braile Biomédica Indústria, Comércio e Representações Ltda.

Av. Pres. Juscelino K. de Oliveira, 1505 • Jardim Tarraf I

Zip Code 15091-450 • São José do Rio Preto • SP • Brazil

Phone #55 17 2136 7000 • Fax #55 17 2136 7040

Federal Tax ID: 52.828.936/0001-09

Technician in charge: Vladimir D. A. Ramirez • CRF/SP 9010

ANVISA Reg. nº 10159030024

www.braile.com.br

BRAILE
BIOMÉDICA