

Organic Tubular Arterial Graft Simple Corrugated Tube



Experience
Quality
Performance


 **BRAILE**[®]
BIOMÉDICA

Description:

The Organic Tubular Arterial Graft Simple Corrugated Tube is made with a bovine pericardium molded in a special support, treated with glutaraldehyde and conserved in formaldehyde 4%. This treatment gives it proper characteristics of resistance, flexibility and lack of antigenicity.

Indication:

The Organic Tubular Arterial Graft in Bovine Pericardium Simple Corrugated Tube is indicated for reconstruction of aneurysms or aneurysmatic disease, aortoiliac and aortobifemoral revascularization, among others.

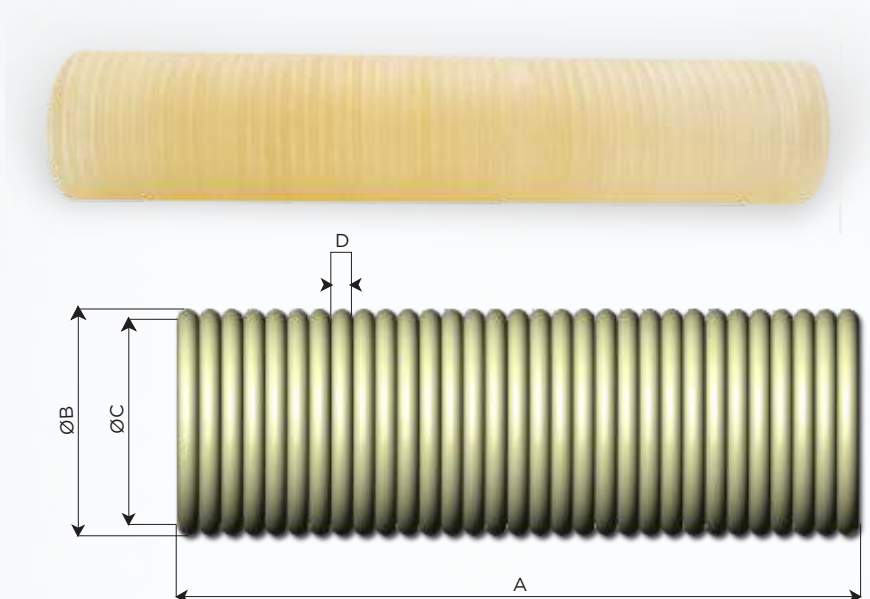

TECHNICAL SPECIFICATIONS (mm)					
Size	A ± 10	ØB ± 0,5	ØC ± 0,5	D ± 0,5	Code
19	120	19	17	1,7	6932
21	120	21	19	1,7	6990
23	120	23	21	1,7	7025
25	120	25	23	1,7	7067
27	120	27	25	1,7	7096
29	120	29	27	1,7	7140
31	120	31	29	1,7	7153
33	120	33	31	1,7	7210

Catalogue Code 600214
Rev. 01 - 07/2018

