

# Rigid Ring for Valve Annuloplasty - Gregori



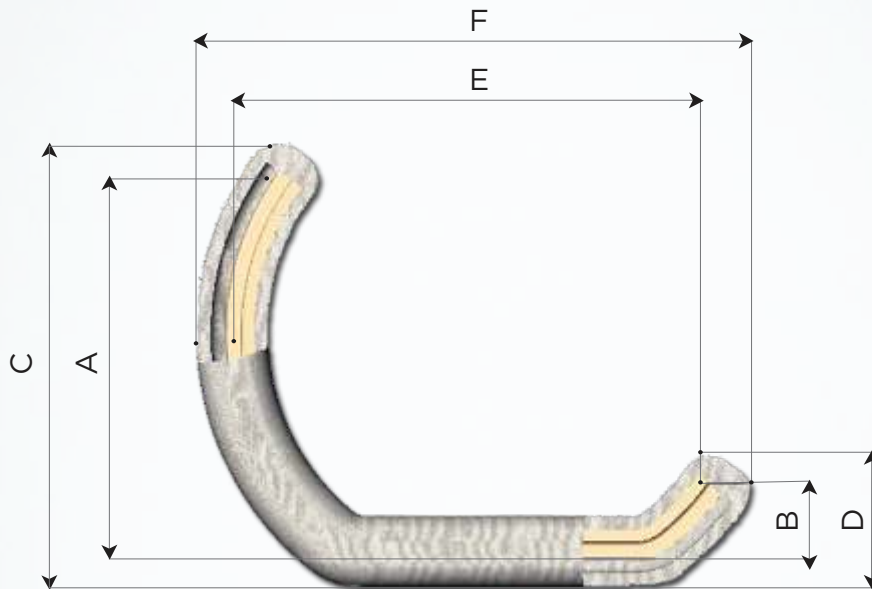
Easy Usability  
Accurate Mitral Ring  
Remodeling Open, permitting  
Mitral Ring Growth  
Experience

**Description:**

A semicircle shape with a rectified portion from its medium point to the next of its end which will correspond, when applying it, to the posterior medium commissure. Made of stainless steel, covered by a silicone tubing and coated by polyester.

**Indication:**

Rigid Ring for Valve Annuloplasty - Gregori is used for Mitral valve repair in reconstructive surgery of Mitral valve, whose purpose is to put the retracted mitral valve in its anatomic natural position, remodel the distended commissures and correct the Mitral ring dilatation, distributing the traction forces by means of wrinkling.

**TECHNICAL SPECIFICATIONS (mm)**

Size	A ± 0,5	B ± 0,5	C ± 0,5	D ± 0,5	E ± 0,5	F ± 0,5	Code
24	17	3	18,5	4,5	20,5	22	2158
26	18	3	19,5	4,5	22,5	24	2161
28	20	3	21,5	4,5	24,5	26	2474
30	21	3	22,5	4,5	26,5	28	2187
32	22	4	23,5	5,5	27,5	29	2190
34	24	4	25,5	5,5	29,5	31	2202
36	25,5	4	27	5,5	30,5	32	2215
38	27	4	28,5	5,5	33,5	35	42147