


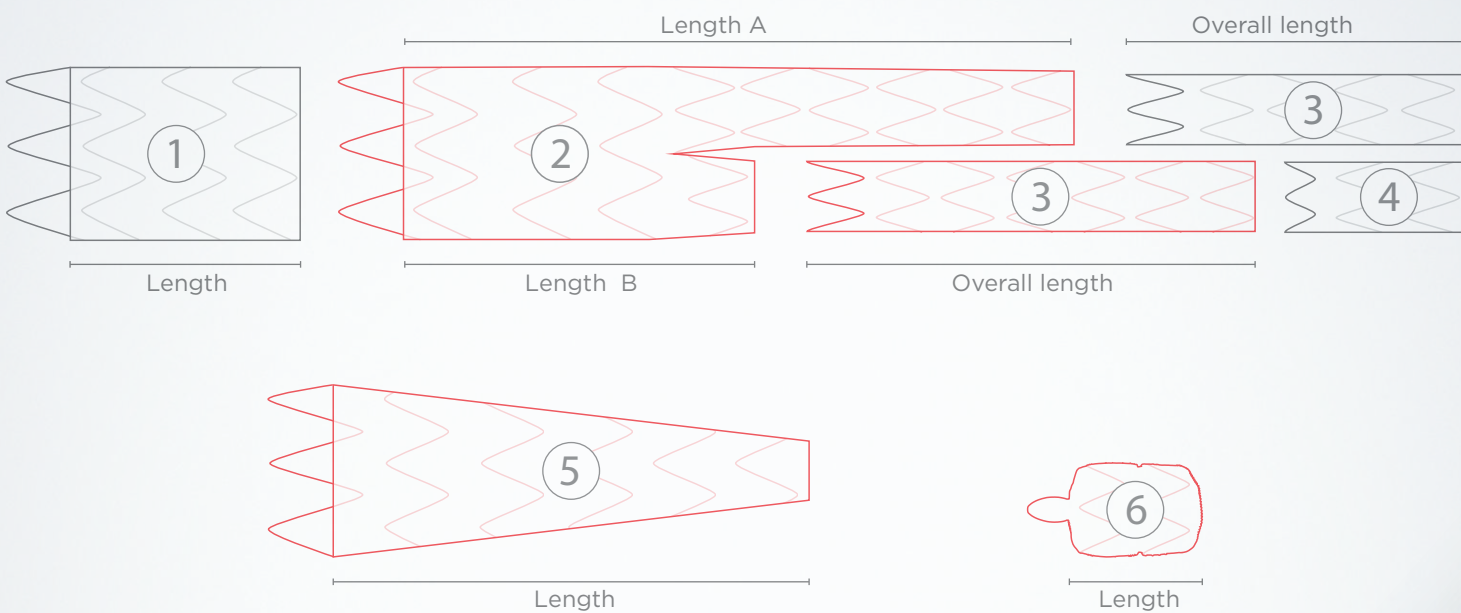

1 Aortic Extension	2 Main Body	3 Iliac Extension	4 Distal Cuff	5 Aortic- iliac	6 Ocluser
Diameter 24 - 36 mm	Proximal Diameter 24 - 36 mm	Proximal Diameter 14 mm	Proximal Diameter 12 a 22 mm	Proximal Diameter 24 a 36 mm	Diameter 14 a 22 mm
Length 40 mm	Distal Diameter 14 mm	Distal Diameter 12 - 14 - 16 - 18 20 - 22 mm	Distal Diameter 12 a 22 mm	Distal Diameter 14 mm	
	Length A 80 - 155 mm	Length 40 - 60 - 80 - 100 120 - 140 mm	Length 40 mm	Length 90 e 150 mm	Length 20 mm
	Length B 75 mm				
Free Flow 20mm	Free Flow 20mm			Free Flow 20mm	
Profile 18 - 20Fr	Profile 18 - 20Fr	Profile 16Fr	Profile 16Fr	Profile 18 e 20Fr	Profile 16Fr

Stent-Graft Dominus
ANVISA Reg. nº 10159030093

***Refer to availability of other dimensions.**

Braile Biomédica Indústria, Comércio e Representações Ltda.
Av. Pres. Juscelino K. de Oliveira, 1505 • Jardim Tarraf I
Zípe Code 15091-450 • São José do Rio Preto • SP • Brazil
Phone #55 17 2136 7000 • Fax #55 17 2136 7040
Federal Tax ID: 52.828.936/0001-09
Technician in charge: Vladimir D. A. Ramirez • CRF/SP 9010
ANVISA Reg. nº 10159030096
www.braile.com.br



LINUS[®]

STENT-GRAFT





Self-expanding metallic stent graft, coated with polyester and Nitinol structure. Attached by Radial Force


► **Free Flow**

Non-coated stent for expansion of attaching area



► **Suture**

Sutures made of polyester for attaching the metallic framework

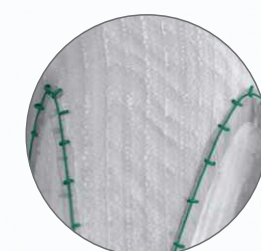

► **Coating**

Low porosity tubular polyester

► **Nitinol**


Superelastic metal alloy with thermal memory
Highly resistant to corrosion and fatigue

► **Ypsilateral**

Exoskeleton

► **Radial Force**

Proximal and distal areas with increased radial force for improving attachment and sealing

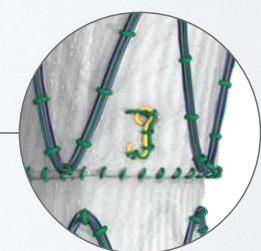

► **Radiopaque Marks**

Proximal and distal areas for easing the location and position of stent-graft

► **Contralateral**

Modular Extensions



► **Balloon Catheter**

Compliant for accommodating stent-graft 12 Fr/45 mm
Supplied separately

Recommended oversize of 15 to 20% concerning the attaching neck diameter
Safety lock