


Tubing Set

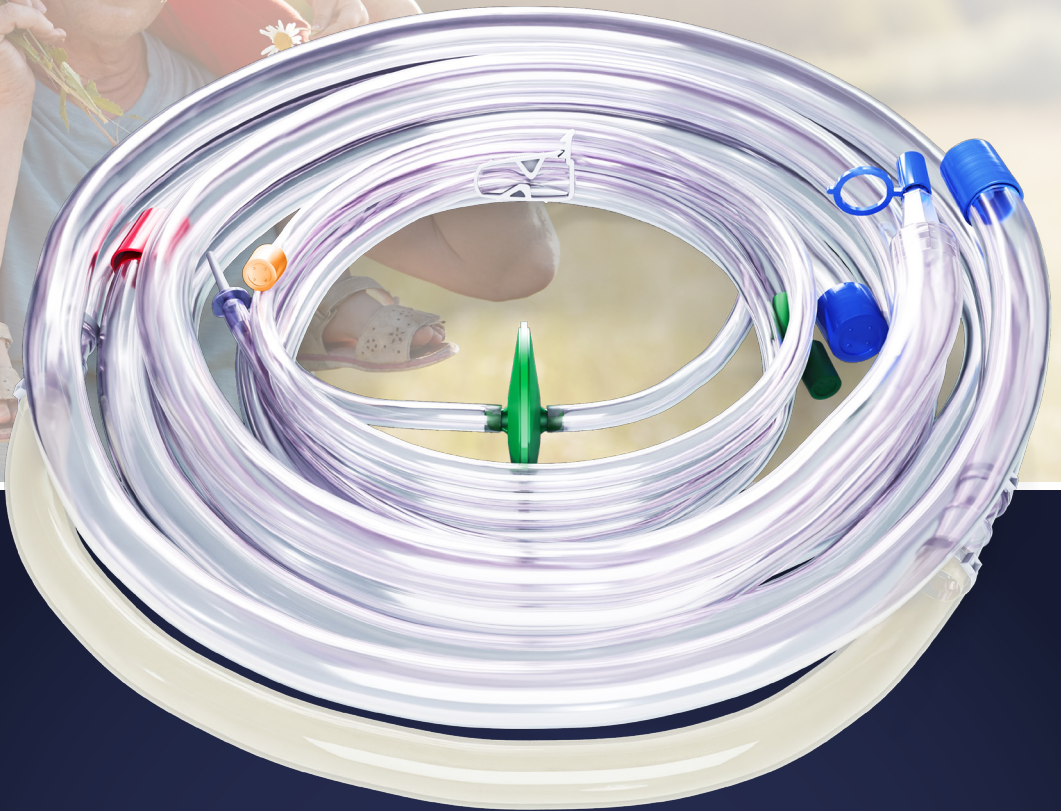

 **BRAILE**
LIFE-SAVING TECHNOLOGIES

CE
2797

The Tubing Set consists of a kit of tubes, connectors and accessories made of medical grade, non-toxic, biocompatible and non-pyrogenic materials. The product consists of the following lines: cava, arterial, pump, aspirators, serum, oxygen, recirculation, arterial filter, extensor and accessory.


Coated models include the biocompatible BRCoating technology on all blood-contact surfaces, forming a biopassive layer designed to reduce blood interaction with the device.

Indication



The Tubing Set is intended to be used for patients who need surgical procedures requiring cardiopulmonary bypass (CPB). The Tubing Set is for use with patients when required blood flow rate (cardiac output) will not exceed:

- Tubing Set Adult: 7 L/min
- Tubing Set Pediatric: 3 L/min
- Tubing Set infant: 2 L/min
- Tubing Set Infant 1500: 1.5 L/min
- Tubing Set Newborn: 0.8 L/min

This device shall only be used by properly trained and qualified personnel.

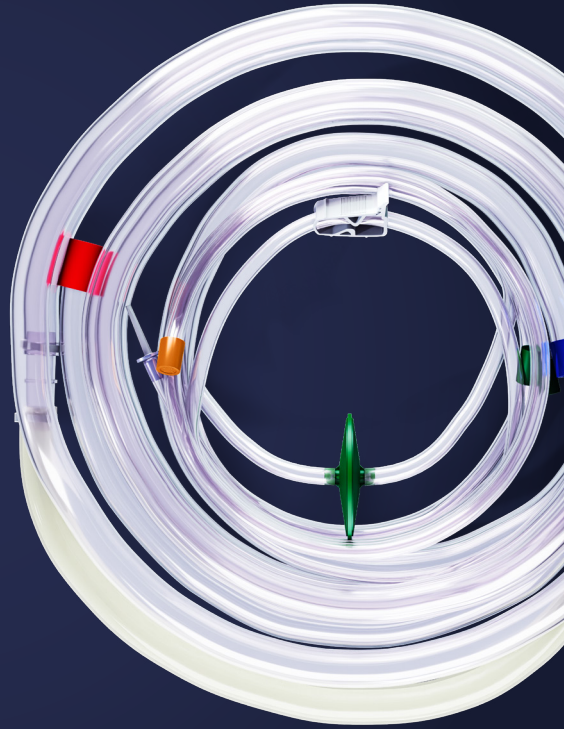


Code	Description	Coating
617322	Tubing Set Adult	No
617323	Tubing Set Infant	No
617324	Tubing Set Infant 1500	No
617325	Tubing Set Pediatric	No
617326	Tubing Set Newborn	No
622286	Tubing Set Adult BRCoating	BRCoating
622290	Tubing Set Infant BRCoating	BRCoating
622294	Tubing Set Infant 1500 BRCoating	BRCoating
622298	Tubing Set Pediatric BRCoating	BRCoating
622302	Tubing Set Newborn BRCoating	BRCoating

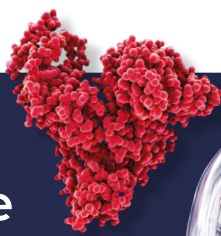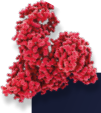

Product covered with biocompatible BRCoating®

A solution made of recombinant albumin which covers the product's inner surface, reducing its interaction with blood.

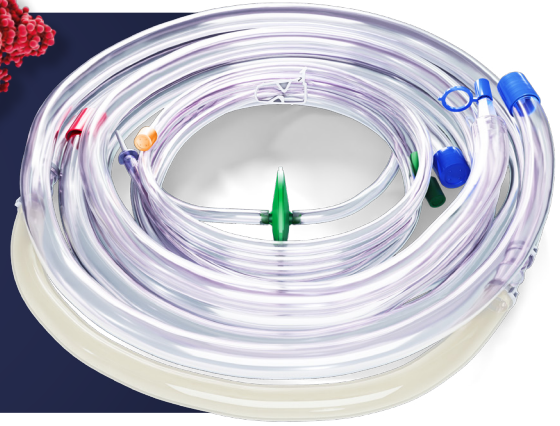

Follow us on social media:

 /brailebiomedica

 @brailebiomedica



Braile Biomédica Indústria, Comércio e Representações Ltda.

Av. Pres. Juscelino K. de Oliveira, 1505 | Jardim Tarraf I

CEP 15091-450 | São José do Rio Preto - SP | Brasil

CNPJ: 52.828.936/0001-09

Technical responsible: Vlademir D. A. Ramirez - CRF-SP 09010