


INOVARE

Alpha



 **BRAILE**

LIFE-SAVING TECHNOLOGIES

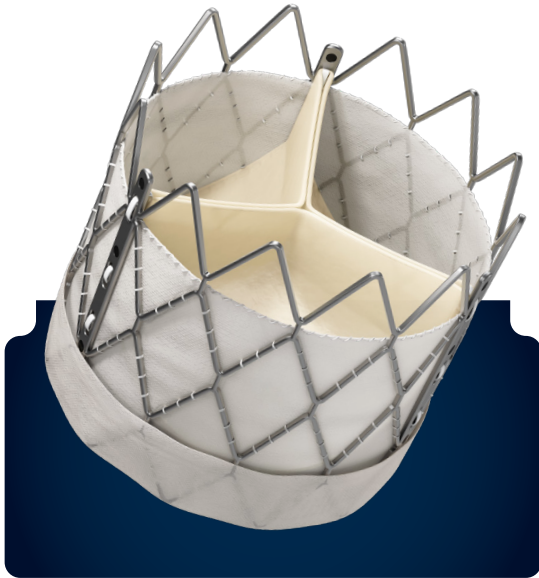

Based on the proven design and performance of the INOVARE valve family.

Indicated for patients with Calcified Aortic Valve Stenosis and for replacing the mitral valve through a transcatheter valve-in-valve implant procedure.

Benefits and Differentiators




Large open cells

Facilitate access to the coronary ostium



Outer polyester (PET) skirt

Efficient sealing to reduce paravalvular leak

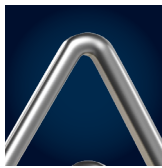

Cobalt-chromium alloy

Structure with increased radial force and reduced recoil




Posts with radiopaque markers

Stability, precision, and safety

- Leaflets made from a single bovine pericardium patch: fewer sutures and larger effective orifice area (EOA).
- Oxidative Treatment – Reox: anti-calcification treatment for bovine pericardium.
- Optimization of surgical time: reduced aortic cross-clamp and cardiopulmonary bypass time.

A solution for patients who are in the gap between conventional surgery and TAVI, patients with a small aortic annulus (risk of mismatch), need for concomitant procedures, and minimally invasive surgeries.

Delivery System

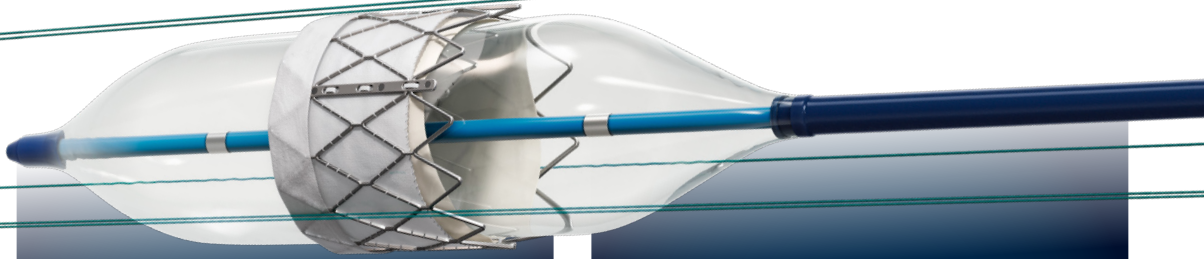


Non-compliant balloon

Optimizes prosthesis positioning in the valvular annulus through uniform expansion

Flexible shaft

Adjustment of angulation and approach to the aortic annulus


More than a valve, a long-term strategy.

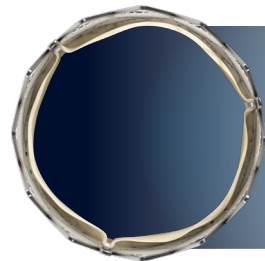


Procedure Components




Inovare Alpha® Biological Valvular Prosthesis
 Sizes 20, 22, 24, 26, 28 and 30 mm

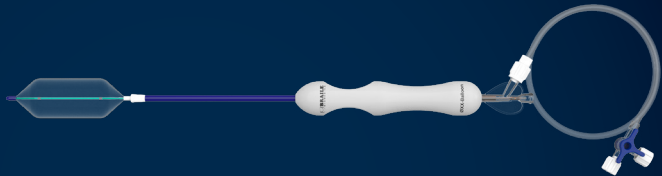

Balloon Delivery Device

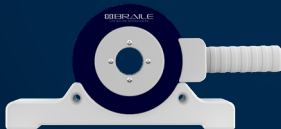

Crimper

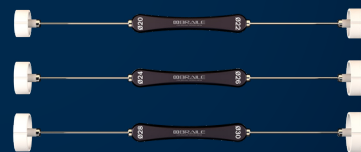

Valve Sizers Kit

Inovare Alpha® Kit

Code	Sizes
604405	20 mm
604406	22 mm
604407	24 mm
604408	26 mm
604409	28 mm
607249	30 mm



For more information, click here or use the QR code beside it.

Follow us on social media:

 /brailebiomedica

 @brailebiomedica



Braile Biomédica Indústria, Comércio e Representações Ltda.

Av. Pres. Juscelino K. de Oliveira, 1505 | Jardim Tarraf I

CEP 15091-450 | São José do Rio Preto - SP | Brasil

CNPJ: 52.828.936/0001-09 | Anvisa no. 10159030095

Technical responsible: Vlademir D. A. Ramirez - CRF-SP 09010